# Gebruiksaanleitung Manual Mode d'emploi

(D)

Seite 2


(GB)

Page 4

(F)

Page 6

**Alpha<sup>®</sup> Antigleitkissen**  
**Alpha<sup>®</sup> Non-Slip Cushion**  
**Alpha<sup>®</sup> Coussin Anti-Glissement**

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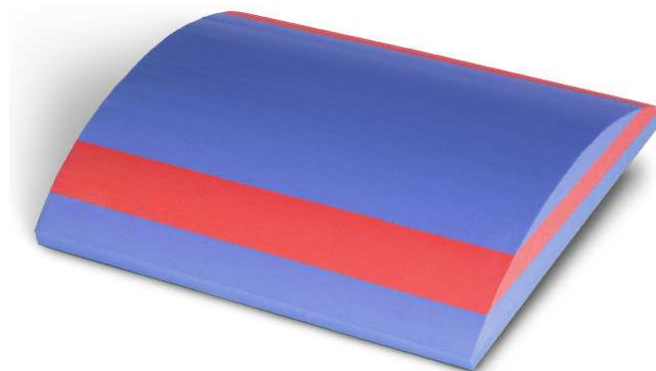## 1. Specifications and Scope of Supply

<b>Item no.:</b>	<b>PM-8050</b>
<b>Dimensions:</b>	50 x 40 cm
<b>Material:</b>	Polyethylene, latex free
<b>Colour:</b>	blue/red
<b>ILN:</b>	4250355380500
<b>Storage:</b>	-5 – 40° C protect from sunlight

**Scope of supply:**

- Non-slip cushion
- Manual

**Made in Germany**


max.  
200 kg

CE
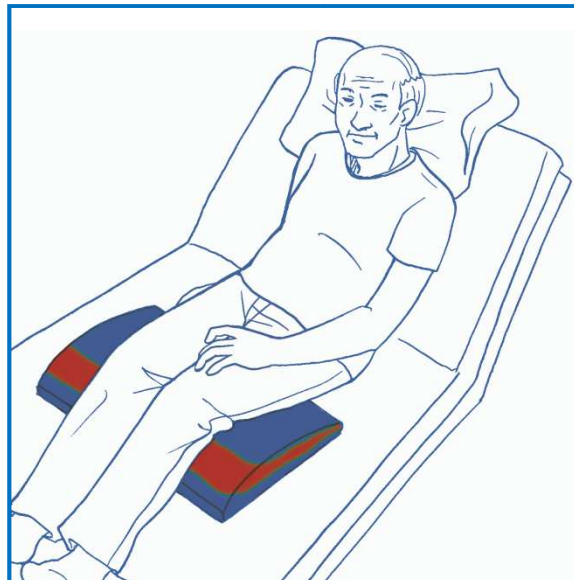
## 2. Use / Placement

The non-slip cushion reduces the risk of a patient seated in bed slipping down. It is placed underneath the patient's thighs.

**This product is suitable for re-use.**

## 3. Maintenance / Care

The product does not require any special maintenance. You can use a normal household cleaner and wipe it down with a moist cloth regularly for reasons of hygiene. A wipe-down disinfection can be performed with the usual disinfectant solutions and isopropyl alcohol. There is also some care information on the label.

## 4. Warranty

- Petermann gives a warranty that the product is devoid of material defects and manufacturing defects.
- The warranty is valid for a period of **3 years**. It does not cover any defects which are due to improper use.
- Petermann cannot be held liable for any accidental or indirect injury or damage.

## 5. Visual and Function Inspection / Safety Information / Re-use

---

Repositioning products are also subject to natural aging. Check the equipment regularly and carefully to make sure it is working properly and especially before re-use. Make sure that the product is in good condition and that the material is not damaged. The high-friction has a rough surface to prevent slippage. The best way to test these properties is by comparing the product with a new product. For re-using the product it must be disinfected before, see 3.

### **WARNING:**

- ***Faulty products must no longer be used.***

## 6. Product Labeling

---

The Petermann Alpha® non-slip cushion is a medical product, class I and complies with the MDD directive 93/42/EEC.

Petermann GmbH is certified in accordance with the directive DIN EN ISO 13485 by the German TUV.

## 7. Disposal Information

---

For the non-slip cushion there are no special disposal rules; it can be disposed of in accordance with national and local regulations as necessary.

If you have any questions about the use of our products, please call our service number:

**+49 (0) 9868 9339-0**

Petermann GmbH • Ziegelhaus 12 • D-91601 Dombühl  
Postal address: Petermann GmbH • Waldweg 3 • D-91601 Dombühl  
Phone: +49 (0) 9868 9339-0 • Fax: +49 (0) 9868 9339-20  
E-mail: petermann@pm-med.de • www.pm-med.de

"Information subject to technical modifications to our products and deviations in colour or printing errors"

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Petermann GmbH est certifié par le TUV conforme à la directive DIN EN ISO 13485.
Adresse postale: Petermann GmbH • Waldweg 3 • D-91601 Dombühl